

 KlarityATHENA™

SUPINE BREAST SYSTEM



KlarityATHENA™

BREASTBOARD + WINGBOARD

The Klarity Athena Supine Breast System offers flexible solutions for better customized positioning.

- Positioning options include:
- Vacuum Bags
 - Klarity Cushions™
 - S-type Head Masks
 - Head & Thorax Masks

Supports up to 440 lbs

Elevations at 5°, 10°, 15° and 20°



R610-PDCF
Carbon Fiber
Athena Breastboard

Arm Supports and Hand-grips not included with the R610-PDCF

KlarityATHENA™


BREASTBOARD

The Athena Breastboard can also be equipped with adjustable arm cups and hand grips as an alternative to the Wingboard.



R610-NDCF
Carbon Fiber
Athena Breastboard
shown with Arm Support

The Athena Supine Breast System includes a special headrest that can facilitate both medial and lateral head placement.

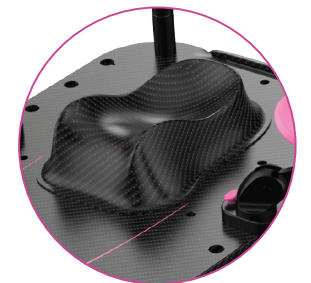

KlarityATHENA™

WINGBOARD

The Athena Wingboard is also available as a separate baseplate. It attaches easily to the Breastboard and can be used as a standalone device.



R632-NTCF
Carbon Fiber
Athena Wingboard

KlarityATHENA™


SUPINE BREAST SYSTEM



The Athena is made of lightweight carbon fiber. MR Compatible versions are also available.



Both the Athena baseplate and the Athena wingboard are compatible with all S-type head masks.

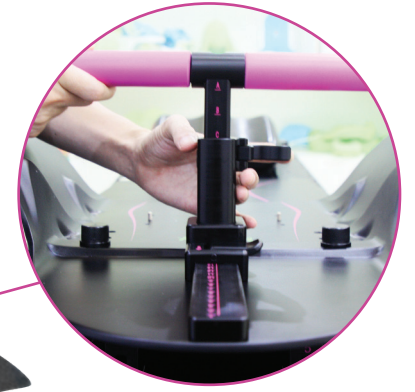

The T-handle can be adjusted both horizontally and vertically at multiple heights.



Elevations at 5°, 10°, 15°, 20°.



Adjustable hand grips offer comfortable and reproducible lower arm placement for one-arm overhead positioning.



The V-shape at the hips relieves pressure to improve both comfort and stability. A 70cm bore can be accommodated when elevated at 5°.


The Klarity Athena Supine Breast System offers flexible solutions for better customized positioning.



Klarity Athena™ Supine Breast System
+ S-Type Head Mask
+ Vacuum Bag



Klarity Athena™ Supine Breast System
+ Head and Thorax Mask (open at treated region)
+ Vacuum Bag

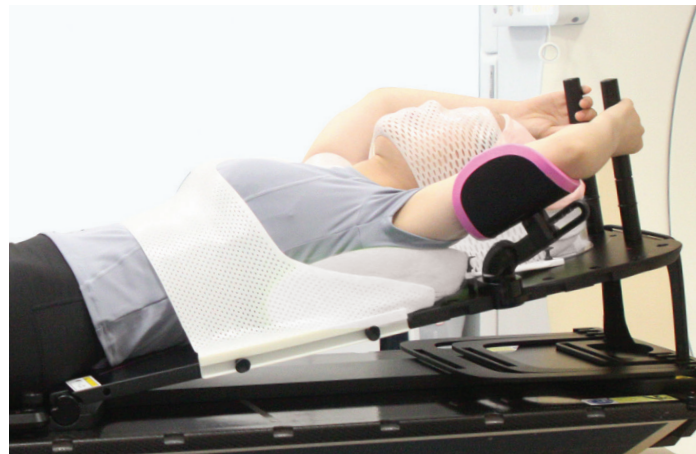

Klarity Athena™ Breastboard
+ Head & Thorax Mask
+ Klarity Cushion™
+ Two-arm overhead positioning with upper hand grips



Klarity Athena™ Breastboard
+ Klarity Cushion™
+ One-arm overhead positioning with upper and lower hand-grips


COMPATIBLE ACCESSORIES



Vacuum Bag:
R7661-20C Vacuum Bag for Upper Body Support
58.5x90 cm
20 Liter Fill



Klarity Cushion™
R550-U3
23 x 35 cm

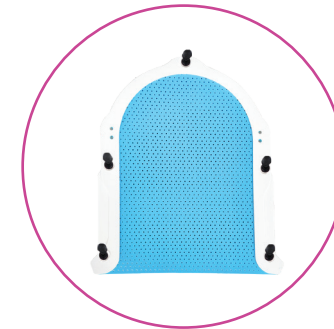

S-type Head Mask
Blue S-Type Head Mask
10% perforation
3.2mm



S-type Head Mask:
Green S-Type Head Mask
36% perforation
2.4mm

Compatible with all standard S-type masks.
Please see our website for a complete selection.

ADDITIONAL ACCESSORIES



Vacuum Bag:
MR7677-25N



Head & Thorax Mask:
R322-S424C1-L (White, Left side)
R322-S424C1-R (White, Right side)

THE ATHENA IS ALSO AVAILABLE AS MR COMPATIBLE



MR R610-PDKF
MR Compatible Athena Breastboard and Wingboard



MR R632-NTKF
MR Compatible Athena Wingboard



Klarity Medical Products USA
600 Industrial Parkway, Heath, Ohio 43056

Toll-free: 877-599-8905 Tel: 740-788-8107 Fax: 740-788-8109

sales@klaritymedical.com www.klaritymedical.com